



Eye Care Instrument

User's Manual



Born for love • Living for Beauty

Welcome to use inFace

inFace eye care instrument is independently developed and manufactured by Shenzhen inFace Technology Co. Ltd., which is affiliated to Shenzhen Migoo Smart Co. Ltd. Since the company was founded, it has been designing the core concepts of practicality, functionality and safety, aiming to create products with technology and intellectual, and striving to explore and innovate.

inFace eye care instrument

Rechargeable, negative ion, vibration and hot compress eye's care.

Function of the product

- Lifting and tightening · Improving fine lines · Diminishing eye circles
- Smoothing eye bags · Nutrient introduction · No allergens

Product Features

Mode 1 white LED is vibration massage: relieves skin pressure and promotes metabolism.

Mode 2 red LED is warm and vibration: warm pores, introduce nutrient essence, and promote absorption.

Mode 3 green LED is vibration plus negative ion: negative ions interact, nutrients directly reach the bottom of the muscle.

Precautions


Please read the user manual carefully before using the product and keep it for future reference.

1. Do not change the use method indicated in this manual without authorization.
2. The instrument is prohibited in the following areas: eyeball, upper eyelid, larynx, wound, skin disease, acne inflammation, rash, itchy redness, etc.
3. Pregnant women are not allowed to use the ion function, please follow the doctor's advice for other functions and use with caution. Not to use the instrument during menstrual period and on anaphylaxis skin.
4. If you wear contact lenses, please remove it before use.
5. Before using the instrument, please clean the skin that needs to be treated. If there is makeup, please use it after removing makeup.
6. When using, the face should be fully moisturized. If the skin becomes dry and difficult to slide the instrument, please add skin care products. Do not slide barely under dry conditions.
7. Relax, hoarse, spotted skin needs a fixed massage, do not pull the skin
8. When using this instrument, do not stay at a certain surface of the skin all the time, but need to move gently.
9. When using the instrument for the first time, please observe the skin condition at any time.
10. This instrument is not waterproof. Do not use the instrument in a bathroom or other humid place, or charge it, or place it to avoid damage caused by water ingress.
11. If the instrument makes abnormal noise, stop using it immediately and contact the sale or service point for repair (please bring the warranty card).
12. Do not use alcohol to clean the instrument.
13. Attention must be paid to the use of heat-insensitive people.

Accessories and appearance

Open the box and you will see:

1. inFace eye care instrument x 1
2. DC charging cable x 1
3. Storage bag x 1
4. User's manual x 1

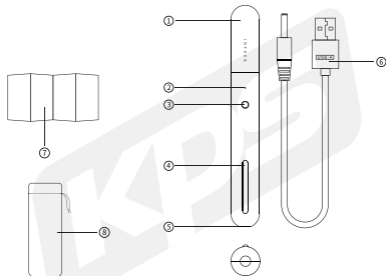

- ① Massage head
- ② LED indicator
- ③ Switch button
- ④ Ion sensing strip
(touch for the ion function)

- ⑤ DC charging interface
- ⑥ DC charging cable
- ⑦ User's manual
- ⑧ Storage bag

Product pavameters

Product Name	inFace eye care instrument
Specification Model	MS5000
Product dimension	131.4mm*18.1mm*17mm
Product net weight	about 46g
Working time (tested condition is DC3.7V)	Mode 1 vibration function is about 250 minutes of continuous work.
	Mode 2 vibration + warm function for about 30 minutes of continuous work.
	Mode 3 vibration + ion function is about 200 minutes of continuous operation.
Warm mode	$39 \pm 1^{\circ}\text{C}$
Battery capacity	3.7V/500mAh
Charging current	250mA
Charging voltage	USB DC5V
Charging time	about 2 hours
Standby current	$<30\mu\text{A}$
Working temperature	$0\sim 40^{\circ}\text{C}$
Material composition	Copper, aluminum alloy, stainless steel

Note: The figures and related product features described in this manual, and some definitions of the terms, can only be used as a reference. Please be subject to the actual products. In view of the consistent policy of continuous improvement of products, the company reserve the rights to change the product specifications without prior notice.

Method of operation

1. Long press 1.5 seconds to switch on and off;
2. The Product under operation, Double-click the button to switch normal mode and detection mode;
3. The Product under operation, short click the button to switch three functions: Mode 1 (white LED light), Mode 2 (constant temperature $38.5 \pm 1.5^{\circ}\text{C}$, red LED light), Mode 3 (green LED light)
4. The product memorize last operation mode, when switch on will operation with last mode function;
5. In the detection mode, the finger needs to touch the ion conductor stripe, and then the massage head contacts the skin to have conductive loop to work;
6. When not switch any mode for 180 seconds, it will shut down automatically.



Long press the switch



Three adjustable levels



Ion circuit

Tips for use



Take a proper amount of eye cream with the ring finger, apply it from the inner corner of the eye to the temple, turn on the product (It is recommended to use Mode 3 green light).



Move the massage head from the inner corner of the eye along the lower eyelid to the temple.



Do the point-like press as the figure and stay at each point for about ten seconds, repeating 2-3 times.



Follow the direction of the arrow in the figure to continue to circle around the eye, usually 3-5 turns.



Apply a proper amount of skin care products or essence to your face and turn on the product according to your needs. (It is recommended to use Mode 1 white light).



Touch the ion conductor stripe and move the massage head slowly from the face along the arrow directions in the figure.

Charging instructions

1. This instrument is USB in-line charging, charging time is about 2 hours, for the working time of each function, please refer to the product parameters.
2. When the battery capacity goes low, the orange light is in a state of respiratory flashing. When the battery is fully charged, the green light is always on, and the instrument cannot work during charging.
3. The instrument will make " di di di " sound when the battery is low to remind users.
4. The charging cable supplied with this instrument can only be used for charging this product. Do not use it for charging other devices.
5. When the charging is abnormal, please turn off the power immediately.
6. This instrument can only be charged with 5V/1A DC adaptor.
7. Please use and store the charging cable properly, do not fold, stretch, twist, entangle, etc.

Preservation method

1. If it is not used for a long time, please wipe the product clean, put it into the storage bag and store it in a dry and ventilated place.
2. The instrument should be protected from direct sunlight and no corrosive chemicals around it.
3. Keep out of reach of children.
4. Do not place the instrument below zero or above 40 degrees to avoid affecting battery life and accidents.

Common problem

If you find any problems during use, please refer to the following table to check the abnormality. If the problem cannot be solved, please contact after-sales service.

Phenomenons	Solutions
Automatic shutdown after using for 3 minutes	This product has a 3-minute intelligent timing protection function, and the instrument will automatically shut down without any operation for 3 consecutive minutes.
No vibration when you use	Please check whether the following conditions exist or not: 1. If the skin is dry, the instrument can not detect the moisture of the skin, there will be no function, please replenish the skin before using 2. Please check whether the battery is fully charged or not, and whether the LED light is normal 3. When the instrument is operating, the ion sensing strip needs to be held by hand, and the guide head is close to the skin to start the vibration operation.
Redness after using the instrument	There is a micro-current that activates cells through the skin. It is normal for the skin to become red after use.
Can not turn on	The instrument can be operated by switching on and off for 1.5 seconds. If there is no response after 1.5 seconds of long press, please check if the following conditions exist: 1. Is the power sufficient? Please re-operate after fully charged. 2. Whether the LED light works during charging.
Why is it tingling on the skin when using the instrument?	When using, the product will release micro-current, and if it feels tingling, it means that it has a micro-current effect and is a normal reaction.
No response when charging	You can follow the steps below to check it: 1. Check and exclude if the power supply is normal. 2. Is the power supply wire and connector damaged? 3. Does the charging LED have a blinking prompt?

About warranty

1. The product is guaranteed for one year from the date of purchase.
2. For products that are out of warranty, the company provides maintenance services, and customers are required to pay maintenance fee and material costs.
3. If you find that the product has quality issues, please contact the after-sales customer service within seven days from the date of purchase.

Warranty restrictions

1. Failure to use or maintain the product in accordance with the product instruction manual.
2. This instrument is made of metal. When storing and using, please try to avoid the friction, drop or strong external impact of the hard object, so as to avoid the appearance damage, the function damage caused by human damage, is not protected by the warranty.
3. Charging cables, storage bags, storage boxes and other accessories and gifts are not included in warranty service.
4. Damage caused by force majeure such as natural disasters.
5. Product failure or damage caused by human factors (such as incorrect operation, strong impact, scratches, disassembly, bumps, use of substandard power supplies, etc.).
6. This instrument contains precision technology and it is forbidden to disassemble and repair the instrument by yourself. The company is not responsible for the irreparable damage caused by the destruction of the instrument or human damage, and the company will not be responsible for the obligations stated in the after-sales warranty.



WARRANTY CARD

MODEL NO. _____
PURCHASE DATE _____
NAME _____
ADDRESS _____
TEL _____
PRODUCT NAME _____
DISTRIBUTOR _____
ORDERR NUMBER _____
DEFECTIVE PHENOMENON _____

(Fill in the warranty card details when applying for the repairing service.
The warranty period is one year for the whole product.)



Product Name: inFace eye care instrument
Product Model: MS5000
Production Date: _____
Quality Inspection: _____

Shenzhen Migoosmart Co., Ltd.





Please Scan the QR code to
re-view the User Manual

After sales contact

Phone: 0086-400-007-8388

Service Wechat : migoosmart

Working hours: 9:00-18:00, Monday to Saturday

(China time, Day off in holidays and weekend)

Shenzhen Migoosmart Co., Ltd.

Address: Rm.807, Building U6, Junxiang U8 Intelligent Manufacturing
Industrial Park, Gushu Community, Xixiang Sub-District, Bao'an
District, Shenzhen, China

Email: service@migoosmart.com

Zip code: 518000

※ The final interpretation of this manual belongs to Shenzhen Migoo
Smart Co., Ltd.



inFace Official Account



inFace service account