


inFace RF Beauty Instrument

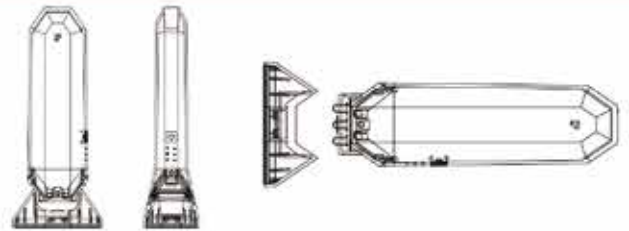

MS6000 RF Beauty Device:

1. Anti-aging
2. Medical_level Titanium Head
3. Performance proved by SGS test result
4. Multi-pole RF interaction

Product Specifications

Product Model:	MS6000
Product Name:	inFace RF Beauty Instrument
Product Size:	154mm*63mm*39mm
Net Weight:	About 125g
Battery Capacity:	1500mAh
Working Voltage:	DC 3.7V
Working Current:	DC3.7V_1000 mA
Rated Power:	8W
Standby Current:	<3uA
Working Frequency:	1M Hz
Charging Mode:	USB DC 5V
Charging Current:	600mA
Charging Time:	About 2.5h
Continuous Working Duration:	<60 Minutes
Working Temperature Range:	0°C~40°C
Material/ingredient:	Titanium Head, ABS, PC, Silicone

inFace RF Beauty Instrument



Operation:

1. Long press the button for 1.5 seconds to turn on with a beep sound.
2. Short press to switch three strength levels, high-medium-low, with beeps and white LED indicator
3. In the on state, long press the button for 1.5 seconds to turn off, with beeps.
4. With intelligent timing function, the instrument can automatically shut down after 10 minutes when it is on.
5. With automatic memory function, the instrument can automatically store the last usage level as the default for the next startup.



Long press the button



Three level adjustable

User Method:



Moving from the collarbone to the mandible (Avoid the throat)

Circling from mandible to the ear



Sliding from inner corner of eye to outer corner of eye, improve the eye end and eyebrow tail (Avoid eyeball and eyelid)

Circling from nasolabial folds to the cheekbones



Cautions:

Important!

1. Do not use the instrument for thyroid or lymphatic areas;
2. Do not use the instrument for throat, eyeballs, ears, upper eyelids, temples;
3. Operate carefully when use the instrument for the areas of thin skin and with blood vessels distributed, such as the corners of the eyes, forehead, neck, mouth, and tibia;
4. people with sensitive skin is recommended to use the low level;
5. Within one hour before using the instrument, do not use flammable liquid skin care products such as alcohol, acetone, perfume, disinfectant, etc.;
6. Except for above 5, other cautions, please, refer to handbook.