


# inFace IPL Hair Removal

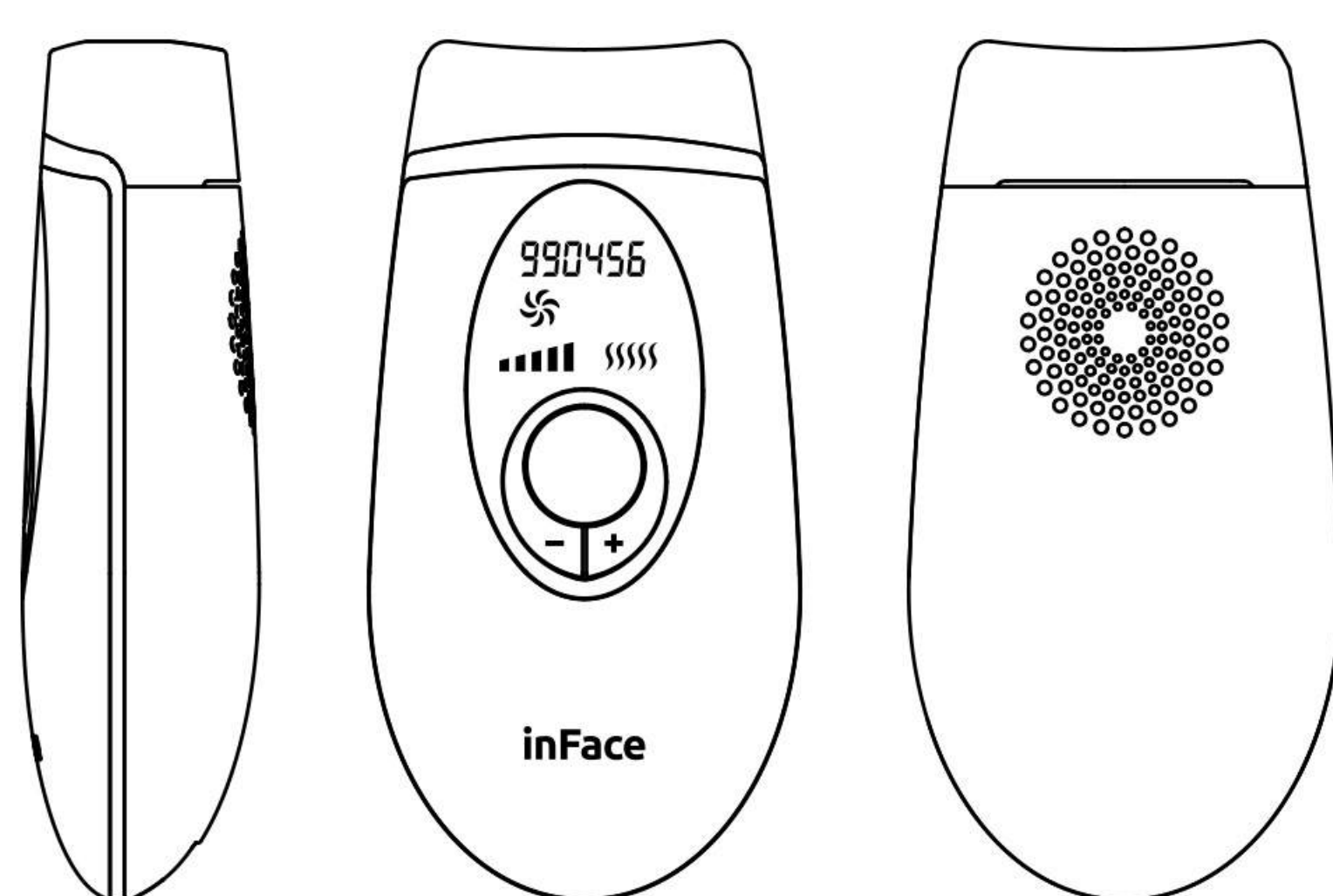

## ZH-01D highlights:

- 1 Million effective flashing , Friendly price
2. Small easy to carry, Large light output area
3. Effective focus , Self-induction light
4. Increase current, Greater the energy

## Product Specifications


Product Model:	ZH-01D
Product Name:	inFace IPL hair removal
Rated Voltage:	12V $\equiv$
Storage Temperature:	0°C~40°C
Operation Temperature:	5°C~35°C
Operation Humidity:	30%-85%
Product Size:	148*75*41mm
Package Size:	240*148*71mm

# inFace IPL Hair Removal



## Preparation:

1. For the parts that need depilation, please use a razor to remove the hair. (Below 1mm)  
(Note: Do not use depilatory cream instead of shaving)

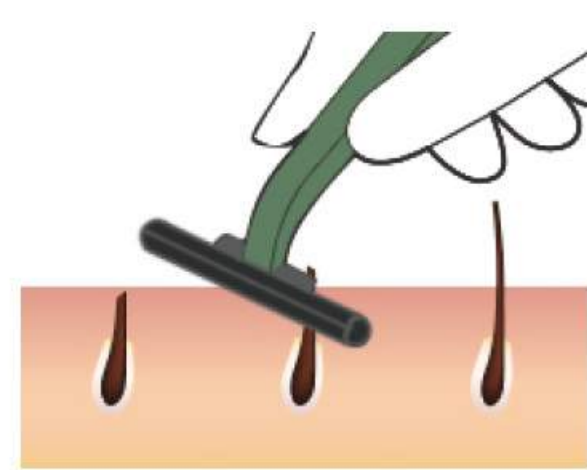

2. Clean the skin to ensure that the skin surface is free of hair, completely dry and free of oily substances.(For example, do not have residues such as fragrance fluid, lotion, perfume, tanning cream and sunscreen)
3. Check whether the device and light port are clean and tidy, and whether there are hair, dust and fibers and other pollutants left. If so, wipe it carefully.
4. Refer to the light intensity table recommended below, adjust the light intensity according to the color of your skin and body hair, and keep yourself in a comfortable state.

### Spectral Reference of Light Intensity Meter

Hair color	White/grey	Red	Beige	Light blond	Brown	Dark brown	Black
Skin color							
White	⊗	⊗	⊗	4/5	4/5	4/5	4/5
Beige	⊗	⊗	⊗	4/5	4/5	4/5	4/5
Light brown	⊗	⊗	⊗	4/5	4/5	4/5	4/5
Brown	⊗	⊗	⊗	1/2/3	1/2/3	1/2/3	1/2/3
Dark brown	⊗	⊗	⊗	⊗	⊗	⊗	⊗
Black	⊗	⊗	⊗	⊗	⊗	⊗	⊗


Note: 1. If you feel any pain or discomfort, please reduce the light intensity;  
2. Recommended light intensity (1-5) In the following table, "X" indicates that this product is not suitable for use.  
5. It is recommended to illuminate depilation parts: legs, armpits, arms, forehead fluff, etc.




## Note:

Male users: please do not use this product on the face with beard or male genital organ.Do not use around eyes, near eyebrows, nostrils, inside ears, or anus.  
Female users: nipples, areola, the labia minora, vagina, do not use.


## Usage Method:

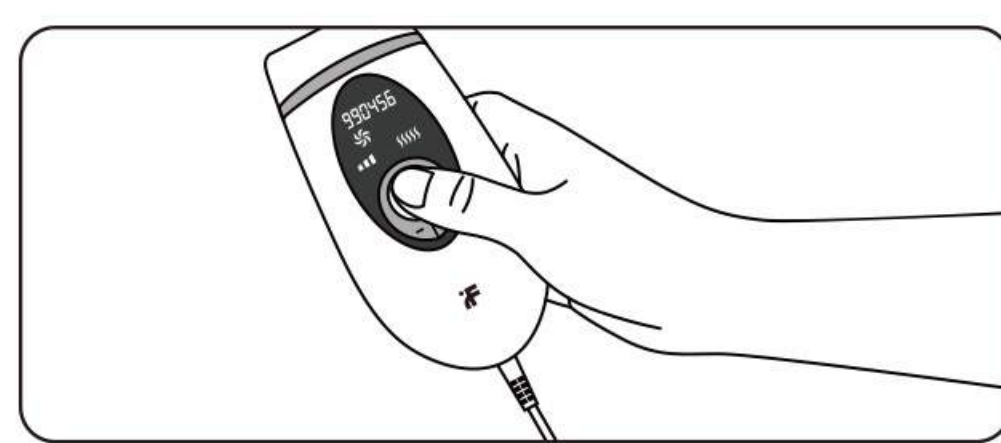


1. Connect the power adapter and long press the key to start the machine.

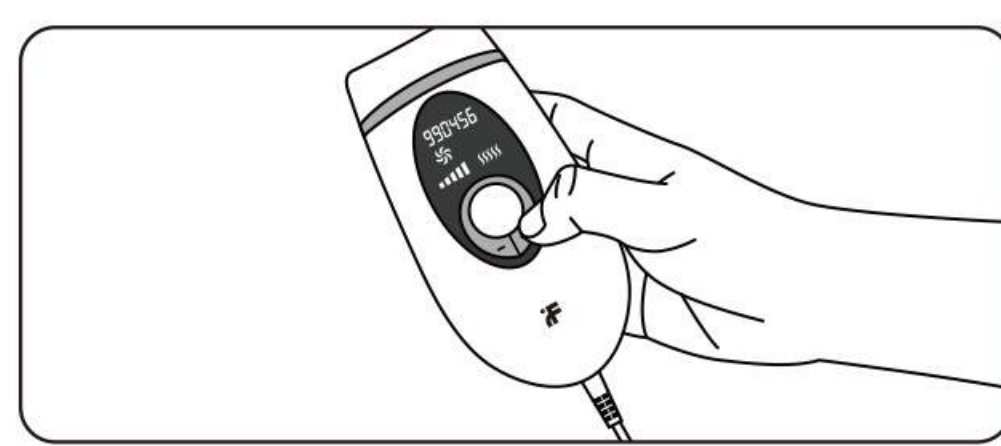


2. Short press the "+" and "-" keys to adjust the light intensity setting; Please choose the suitable level according to your own situation. The higher the level, the better the effect. Novice suggests starting from the lowest level.

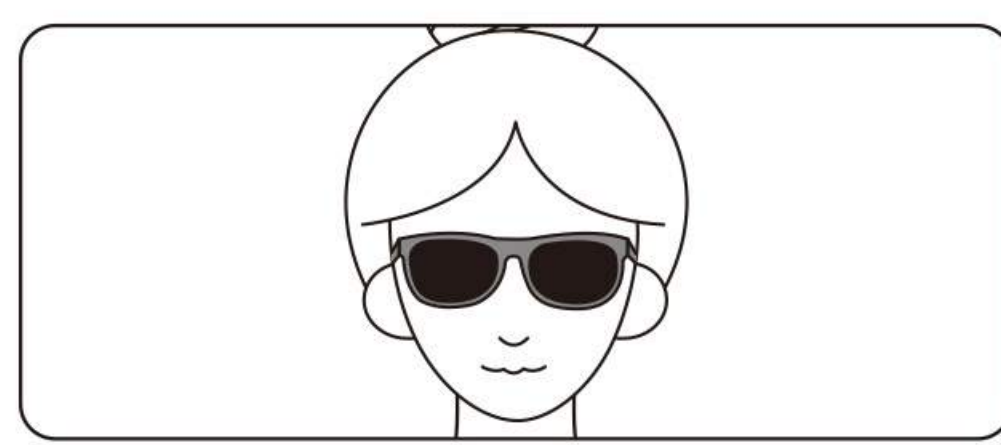

3. When operating, please wear the protective sunglasses which attached in gift box to protect your eyes.

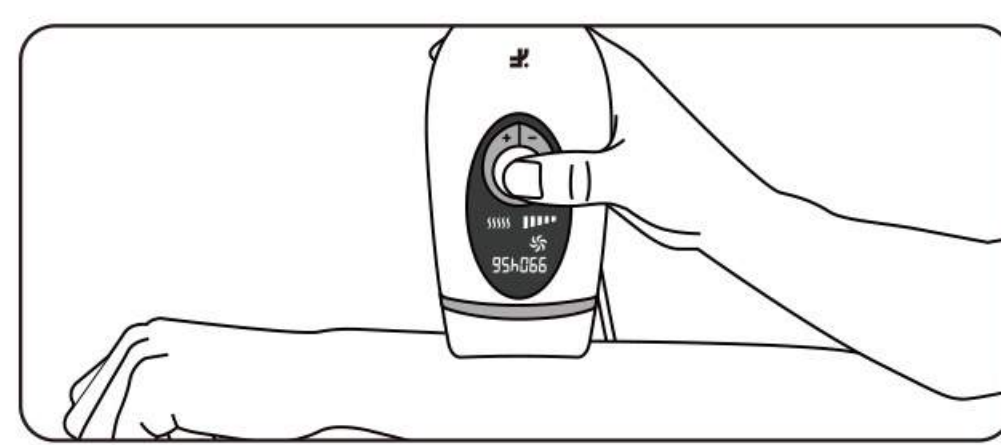


4. Apply this product vertically to the skin, press the key shortly, and begin to illuminate and depilate.

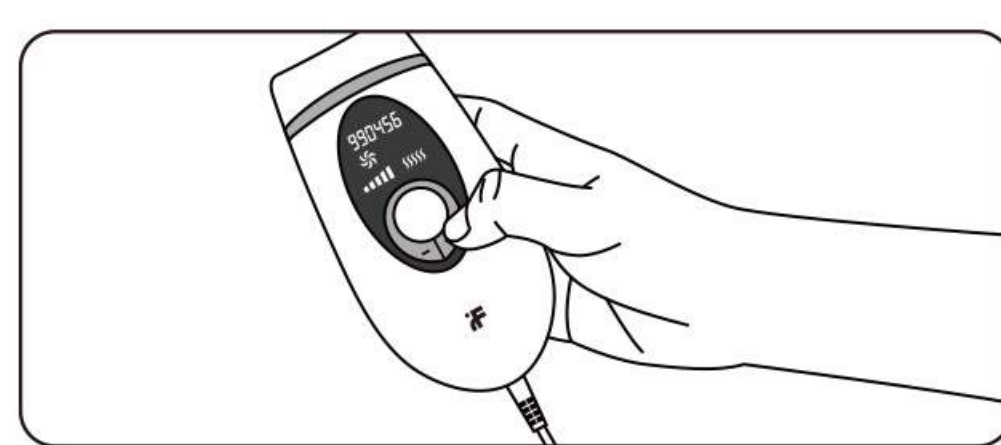

5. The default startup mode of this product is manual mode. Long press the "+" or "-" key to adjust to the

Auto mode (in the automatic mode, this product can continuously illuminate and depilate when contacting the skin, with an interval of 2-3 seconds for each illumination) and move to depilate.

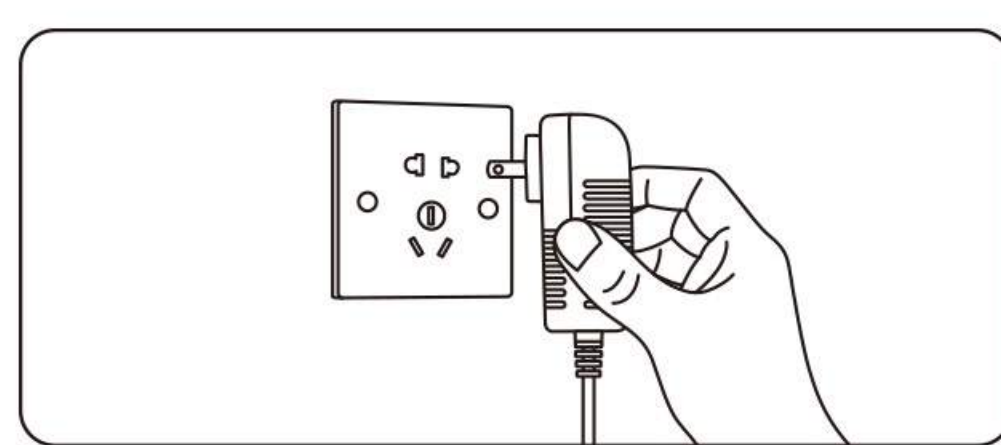

6. After depilation, Long press the key to switch off the product and pull off the power adapter.